
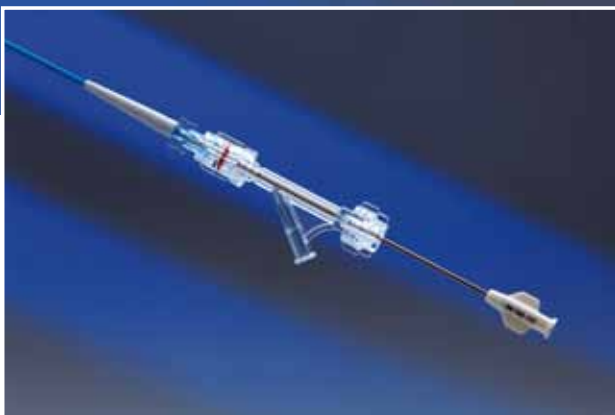


RESISTANT - SFA & Popliteal Stenting

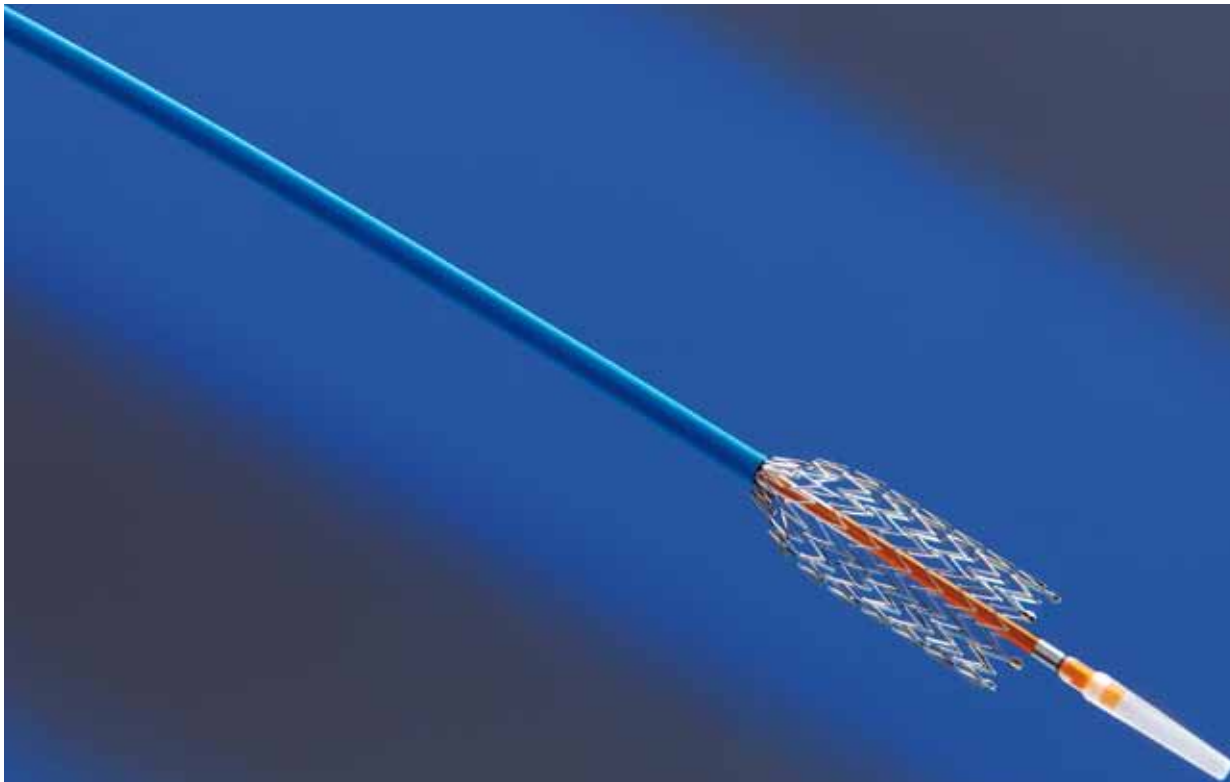

RESISTANT -
is the result of
the optimal balance
of axial flexibility,
bending flexibility
and radial force.

- atraumatic soft tip
- 3 radiopaque Tantalum marker
- 2 additional x-ray marker on shaft
- precise handling and placement of the stent
- approved pull back application device

RESISTANT - SFA & Popliteal Stenting

Stent Details			
Type of design	Slotted tube, 3 interlinks per segment		
Material	Nitinol		
Stent diameter	4-7 mm	8 mm	9-12 mm
Strut thickness	215 µm	215 µm	215 µm
Strut width (main segment)	113 µm	135 µm	160 µm
Strut width (interlink)	111 µm	113 µm	220 µm
Metal coverage	10% - 25%		
Marker	Tantalum		

Stent Delivery System Details	
Usable length	85 cm / 135 cm
Shaft size	6 F
Marker bands	Iridium
Recommended guidewire	0.035"

Order Information								
Stent Diameter mm - Stent Length mm								
		20	30	40	60	80	100	
Shaft- Length 85 cm	5.00 mm	RSO 05/020-085	RSO 05/030-085	RSO 05/040-085	RSO 05/060-085	RSO 05/080-085	RSO 05/100-085	
	6.00 mm	RSO 06/020-085	RSO 06/030-085	RSO 06/040-085	RSO 06/060-085	RSO 06/080-085	RSO 06/100-085	
	7.00 mm	RSO 07/020-085	RSO 07/030-085	RSO 07/040-085	RSO 07/060-085	RSO 07/080-085	RSO 07/100-085	
	8.00 mm	RSO 08/020-085	RSO 08/030-085	RSO 08/040-085	RSO 08/060-085	RSO 08/080-085	RSO 08/100-085	
	9.00 mm	RSO 09/020-085	RSO 09/030-085	RSO 09/040-085	RSO 09/060-085	RSO 09/080-085	-	
	10.00 mm	RSO 10/020-085	RSO 10/030-085	RSO 10/040-085	RSO 10/060-085	RSO 10/080-085	-	
	12.00 mm	RSO 12/020-085	RSO 12/030-085	RSO 12/040-085	RSO 12/060-085	RSO 12/080-085	-	
		20	30	40	60	80	100	
Shaft- Length 135 cm	4.00 mm	RSO 04/020-135	RSO 04/030-135	RSO 04/040-135	RSO 04/060-135	RSO 04/080-135	RSO 04/100-135	
	5.00 mm	RSO 05/020-135	RSO 05/030-135	RSO 05/040-135	RSO 05/060-135	RSO 05/080-135	RSO 05/100-135	
	6.00 mm	RSO 06/020-135	RSO 06/030-135	RSO 06/040-135	RSO 06/060-135	RSO 06/080-135	RSO 06/100-135	
	7.00 mm	RSO 07/020-135	RSO 07/030-135	RSO 07/040-135	RSO 07/060-135	RSO 07/080-135	RSO 07/100-135	
	8.00 mm	RSO 08/020-135	RSO 08/030-135	RSO 08/040-135	RSO 08/060-135	RSO 08/080-135	RSO 08/100-135	

Durability and long term performance RESISTANT



eucatech



competitor A



competitor B

eucatech AG

Gottlieb-Daimler-Strasse 2
D-79618 Rheinfelden, Germany

Tel. +49 7623 71 77 - 0
Fax +49 7623 71 77 - 177

Internet: www.eucatech.de
E-Mail: info@eucatech.de

